

APPLIQUATOR

Fiber-free APPLIQUATOR™
for the precise application of
dental adhesives and liquids.



The Applicator is available
in colour-coded soft and regular
versions and in two sizes (please
see backpage).



transcodent.de

You do precision work.
We create your tools.

FIBER-FREE APPLIQUATOR™

Transcodent has completely re-engineered this patented solution for the application of dental liquids and gels. The Applicator combines a rigid handle with a soft pliable head. The head contains capillary canals that pick up and release precise doses of dental materials. The handle bends and holds its shape for tough-to-reach applications.

Since the Applicator is precision molded, the capillary canals always hold and release the same amount of fluid every time you use one. Instead of painting with a brush, you can apply just the amount of liquid needed for the procedure, reducing waste and the risk of excess material in the restoration site. With a pliable head and no glue sphere, the Applicator flexes into tight spots and corners.

FEATURES

- 100 % fibre - free
- Manufactured as „1 piece“
- Co-molded
- SOFT variant with extra flexible head
- Two sizes: regular / fine
- Designed shape
- Easy to bend (Head holds the position)

BENEFITS

- Precise application
- Bends easily
- Non dripping
- No contamination of the cavity through applicator „fibres“
- No glued parts
- Consistent „transport volume“

FEATURES

- Bonding agent
- Varnish
- Desensitizer
- Cavity liner
- Etchant

PACKAGING


- Trays of 50 pieces
- Trays of 100 pieces



The co-molded head design is optimized to carry a precise dosage of liquid and is available in two different sizes.




Unit Dose

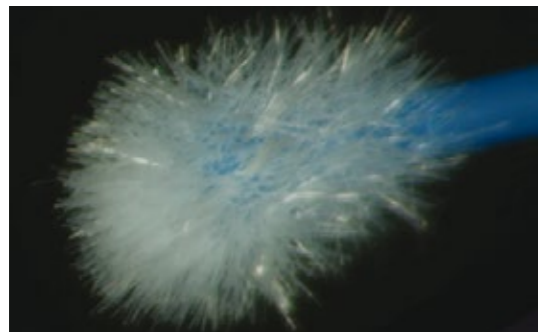

Micro Brushes



FIBERS

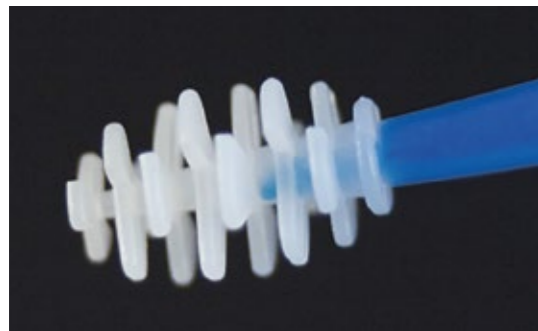

Transcodent APPLIQUATOR



FIBER-FREE



ITEM OVERVIEW

ITEM	SIZE	REF	COLOUR CODE	PACKAGE UNIT
Intro Kit	all sizes	400	all colours	400 (4×100)
regular/standard	regular	401	blue	100
regular/soft	regular	402	green	100
small/standard	small	403	yellow	100
small/soft	small	404	pink	100

OEM variants are available as 50-pcs. and 100-pcs. tray.

The Applicator is a completely new approach to the application of dental liquids and gels. Since it's not a brush and not made from fibers and glue, nothing contaminates the restoration or bonding materials.

The Applicator picks up, holds, and releases the same exact dose every time, giving dentists precise control.



Transcodent GmbH & Co. KG
Zum Kesselort 53
D-24149 Kiel

+49 (0) 431 / 979 979 - 0
info@transcodent.de
transcodent.de

transcodent®